

Breitländer GmbH Mercatorstr. 51 46485 Wesel Deutschland

Tel.: +49 (0)281 319 391 0 Fax: +49 (0)281 319 391 29 E-Mail: mail@breitlander.com

Sample Grinder HK150EXS

Manual pendulum grinder with steel cabinet and integrated spark trap and dust exhauster





HK150 with magnetic plate and special holder for wire samples



Standard vice with medium size jaws



GRINDING MACHINE **HK150EXS**

- GRINDING MACHINE Clean, powerful sample grinding system with stone ready to go
 - Floor standing machine, dustproof, covered working place, halogen spot light
 - Steel cabinet with switch board, dust trap, exhauster and drawer
 - Reproducible flat surface, easily operated pendulum arm in taper bearings
 - Simple and reliable machine construction
 - Specially suitable for hard materials, like white cast irons
 - Adaptable to various qualities by special stone cup wheels
 - Vertical moving of grinding head by handwheel with 0.1 mm scale

Technical Data and Accessory

HK150EXS

www.breitlander.com

Technical Data

Revolutions

Grinding stone


Lighting

Grinder motor 3 x 400 V/50 Hz. 1.1 kW

alt.: 1 x 220 V/60 Hz alt.: 3 x 220 V/60 Hz 2800 rpm / 3440 rpm 12 V / 20 W halogen 150 x 80 x 32 mm

Grinding height, min. 4 mm Grinding height, max. 90 mm

Mass (grinder) approx. 100 kg Mass (cabinet) approx. 50 kg Dimensions (grinder) 400 x 400 x 680 mm Dimensions (cabinet) 600 x 600 x 870 mm

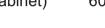

Standard Accessory Vice for sample holding

Grinding stone, tool for stone change Dressing bar for grinding stone Transport handles

Optional Accessory

Magnetic vice with handle and quick exchange arm Special vices - bigger vice for max. diam. 58 mm Hardened vice jaws: 3 mm, 9 mm (regular) and 20 mm Copper cooling blocks

Corundum / silicon carbide cup wheels 150 x 80 x 32 mm Filter bags / exhauster motor protection filter

Spare parts kit: 2 bellows, safety glass, halogen lamp





Cup wheel, open corundum stone



Magnetic plate with quick exchange arm



Vices for different sample diameters



Copper cooling blocks





Exhauster motor protection



Spare parts kit

Pendulum grinder HK150 in combination with disc grinder HK350